

# Norepinephrine ELISA kit

## Ultra Sensitive – Any sample

Ref: BA-E-5200R

This research Norepinephrine ELISA kit enables the ultra-sensitive measurement of Noradrenaline (Norepinephrine) in any biological sample. Intended for research use only, the assay features a lower limit of detection at 1.3pg/ml and a minimal sample volume as low as 1µl.

<b>Sample type</b>	1 - 750µl, all biological fluids
<b>Samples/kit</b>	40 samples in duplicate
<b>Sensitivity</b>	1.3pg/ml (depending on sample volume)
<b>Range</b>	0.2 - 32ng/ml
<b>Assay time</b>	Sample extraction, acylation and enzymatic conversion (4h), ELISA overnight
<b>Reactivity</b>	Reacts with all species
<b>Citations</b>	<a href="#">30+ papers</a>

## INFORMATIONS

---

### Product overview

<b>Product name</b>	Norepinephrine ELISA kit - Research
<b>Description</b>	Enzyme immunoassay (ELISA) for the ultra-sensitive determination of Norepinephrine in biological samples. <b>Research use only</b>
<b>Format</b>	96-well plate
<b>Samples</b>	Any sample
<b>Sample volume</b>	1 - 750µl
<b>Reactivity</b>	Reacts with all species
<b>Standard range</b>	0 / 0.2 - 32ng/ml
<b>Sensitivity</b>	1.3pg/ml
<b>Specificity</b>	No significant cross-reactivity was observed with Noradrenaline analogs such as Dopamine, Adrenaline, Metanephrine, Normetanephrine, Methoxytyramine, 3-Methoxy-4-hydroxyphenylglycol, Tyramine, Phenylalanine, Caffeinic acid, L-Dopa, Homovanillic acid, Tyrosine and 3-Methoxy-4-hydroxymandelic acid
<b>Assay time</b>	Sample preparation 4 hours and ELISA overnight
<b>Storage</b>	Store at 2-8°C for to 6 months
<b>Datasheets</b>	<a href="#">Instructions for use</a> , <a href="#">Material safety datasheet</a>

## PROTOCOLS

---

**Sample preparation**

Sample extraction and acylation (1h45), enzymatic conversion (2h)

---

**ELISA**

Norepinephrine antiserum overnight incubation, revelation and read steps (1h)

---

**Detailed protocol**

[Download instructions for use](#)

## REFERENCES

---

**Product citation**


## Product pictures

---



**Noradrenaline Research ELISA kit**

---

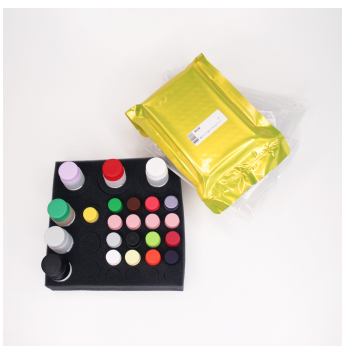

**Noradrenaline Research ELISA kit**

---



**Noradrenaline Research ELISA kit**



**Noradrenaline Research ELISA kit**

## Contact information

Immusmol  
229 Cours de l'Argonne  
33 000 Bordeaux - France  
Tel: +33 (0) 5 6431 1170  
[www.immusmol.com](http://www.immusmol.com)

**To order, review, ask for technical support, visit product page at:**

<https://www.immusmol.com/shop/norepinephrine-elisa-kit-any-sample/>